

BagMixer[®] 400

400 mL lab blenders



interscience

Why choose BagMixer® 400!

Best-seller

- Sold in more than 85 countries
- 35 years experience

High-quality

- Optimal analysis with the sterile **BagFilter®** and **BagPage®**
- In compliance with ISO 7218 standard
- Made in France by **interscience**

Brushless motor

- No wearing parts
- - 50% power consumption
- Extra silent

Safety

- The handle closes the blender and starts the blending
- **Q-Tight®**: no leaks from the top of the bag
- Security drip tray

PATENTED

Click & Clean® removable paddles

- Remove the paddles and easily clean the blending chamber

PATENTED

Adjustable blending power

- Up to 28 kg of pressure: the hardest samples don't resist!
e.g.: nuts, pills, seeds, dog biscuits etc.



before after



Easy cleaning


- All stainless steel
- Door opens at 270°
- Removable **Click & Clean®** paddles

Guarantee

- Shock absorbers and window door **life-time** guarantee
- 3 year warranty (after registration)

Optimal extraction

- From as little as 30 seconds
- Easy control with the digital screen (speed, time, power)

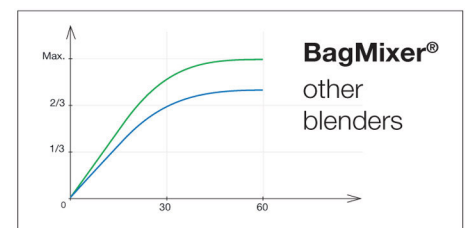

More than 14,000 **BagMixer®** in use every day!

Made in France by **interscience**,
the only manufacturer of blenders and sterile bags for microbiological analysis

Some features are limited to specific models.
For complete technical info, please check online:
www.interscience.com



BagMixer® 400 models:

The **BagMixer®** lab blender guarantees optimal extraction by the action of its two paddles. Adapted to any sample (food, agricultural, chemical, medical, environmental) to be analyzed between 50 mL and 400 mL, it guarantees homogenization in 30 to 60 seconds.



BagMixer® 400 P
 > The essential

BagMixer® 400 W
 > Window door

BagMixer® 400 CC
 > Click & Clean®

BagMixer® 400 S
 > Silent blending

BagMixer® 400 SW
 > Window door
 & LiquidSensor

Ref. 021 230

Ref. 022 230

Ref. 024 230

Ref. 025 000

Ref. 025 100

Also available: packs including bags and accessories; 100 mL and 3500 mL lab blenders. For more info, please contact us.
 For more info on **BagMixer® S & SW** models, please refer to **BagMixer® S & SW** brochure.

BagMixer® packs:

Starter Pack



- 1 **BagMixer®**
- + **Sterile bags:** 500 BagFilter® + 500 BagLight®
- + **Accessories:** 50 BagClip® + 1 BagOpen® + 1 BagRack®

Classic Pack




- 1 **BagMixer®**
- + **Sterile bags:** 500 BagFilter® + 500 BagLight® + 500 BagPage®
- + **Accessories:** 50 BagClip® + 1 BagOpen® + 1 BagRack® + 1 PetriPile® + 1 BagPipet® + 1,000 BagTips® Regular + 6 Clinet® 101

Premium Pack




- 1 **BagMixer®** & 1 **DiluFlow®**
- + **Sterile bags:** 1,000 BagFilter® + 2,000 BagLight® + 1,000 BagPage®
- + **Accessories:** 50 BagClip® + 1 BagOpen® + 2 BagRack® + 2 PetriPile® + 1 BagPipet® + 1,000 BagTips® Regular + 6 Clinet® 101

BagSystem® complete range for sample preparation before analysis

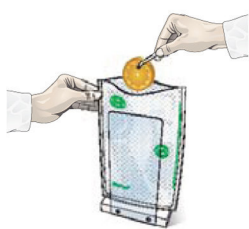

BagFilter®
Collect

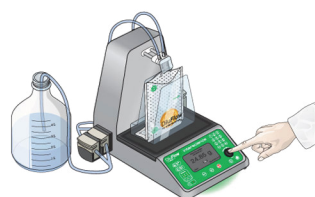

DiluFlow®
Dilute








BagMixer®
Blend



BagTools®
Pipet

Technical specifications

	BagMixer® 400 P	BagMixer® 400 W	BagMixer® 400 CC
Reference	021 230	022 230	024 230
Key features	Removable and autoclavable paddles, Click & Clean® system	-	✓
	Blending power indicator	-	✓
	Adjustable blending power (up to 28 kg of pressure)	-	✓
	Digital screen	-	✓
	Leak free Q-tight® system	-	✓
	Side-by-side paddle stop	-	✓
	Silent brushless motor	-	✓
	Window door	-	✓
	Security drip tray	in option	✓
	270° opening door	✓	✓
	All stainless steel	✓	✓
	Compact and ergonomic	✓	✓
Technical specifications	Door weight	1050 g	950 g
	Blending time	30, 60, 90, 120, 150, 180, 210 sec., ∞	
	Blending speed	fixed speed: 8 strokes/sec.	variable speed: 4, 6, 8, 10 strokes/sec.
	Dimensions (w x d x h)	27 x 38.5 x 28 cm	
	Weight	16.6 kg	
	Power	110-127/220-240 V~ 50-60Hz	
	Capacity	50-400 mL	
Warranties	Made in France	✓	✓
	Life-time warranty	Unbreakable shock absorbers and window door	
	3 year warranty (after registration)	✓	✓
Manufactured in compliance with	    		

Discover our complete range for microbiology



BagFilter®
Sterile bags




DiluFlow®
Gravimetric dilutors



FlexiPump®
Peristaltic pumps



easySpiral Dilute®
Automatic diluter and platers

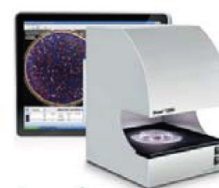

Scan®
Automatic colony counters

Your local distributor

interscience INTERNATIONAL

30, chemin du Bois des Arpents - 78860 St Nom - FRANCE
T: +33 (0)1 34 62 62 61 - F: +33 (0)1 34 62 43 03 - Email: info@interscience.com

interscience USA & CANADA

32 Cummings Park - Woburn, MA 01801 - USA
T: +1 781 937 0007 - F: +1 781 937 017 - Email: sales.usa@interscience.com

interscience CHINA

上海市徐汇区吴兴路277号锦都大厦588室 - 200030
电话: +86 (0)21-64739390 - +86 189 3097 0733 - 邮址: sales.china@interscience.com

interscience SOUTH-EAST ASIA

541 Orchard road - 09-01 Liat Towers - SINGAPORE 238881
T: +65 6933 1389 - +65 8549 1217 - Email: sales.asia@interscience.com

www.interscience.com

Bioswisstec AG Ebnatstrasse 65 8200 Schaffhausen Tel. +41 52 620 33 44