

## Bioswisstec lime green line

- #53250 Bioswisstec 50ml Skirted Conical Tubes, PP, gamma sterilised SAL10<sup>-6</sup>, USP-6 raw material.  
CHF 59.- for 500 pcs.
- #53050 Bioswisstec 50ml Conical Tubes PP, gamma sterilised SAL10<sup>-6</sup>, USP-6 raw material.  
CHF 56.- for 500 pcs.
- #53015 Centrifuge Tubes 15ml PP, with HDPE cap, conical, graduated, with white labeling field (incl. free rack\*), RNA/DNA RNase/DNase free  
CHF 56.50 for 500 pcs.
- #60115 Micro Tube 1,5ml with flat cap, steril, transparent, clear Premium version  
CHF 150.- for 5000 pcs.
- #900025 Bioswisstec Tissue culture flask PS, 25cm<sup>2</sup> Filter cap 0.2µm PE angled neck, working vol. 7ml (also available as 75cm<sup>2</sup> and 175cm<sup>2</sup>)  
CHF 75.20 for 200 pcs
- #800020 Bioswisstec Premium Cell Scraper PP, 22.5cm, turning head, 1.3cm soft PE blade  
CHF 127.- for 100pcs.
- #505050 Cell-Chip, Neubauer "improved" counting chamber, single use with 2 counting chambers, plastic with quartz glass properties,  
CHF 130.- for 50 pcs .




# reliable, compact, high-performance

Take your lab into the **Future**

**HERMLE**  
LABORTECHNIK

### Did you know?

Bioswisstec offers a maintenance service for all devices and machines, so you can focus on your research without worrying about equipment - ask us about options

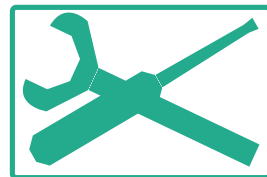

Bioswisstec AG  
Ebnatstrasse 65  
8200 Schaffhausen

bioswisstec.com  
info@bioswisstec.com  
Tel +41 52 620 33 44



premium partner



# The new Centrifuges

## Easy Handling:

The Ergonomic Touch Operation Function enables a smooth and swift operation process.

There is no need to remove surgical gloves, during operation

Strictly a HERMLE application, you can store up to 99 different programs.

For fast centrifugation, there is a quick bottom option. The centrifuge spins, as long as you press the button.

For re-assurance that your samples are finished, an acoustic signal will sound when the centrifugation is complete (This function can be de-activated!).

## Safety

The Z 207 A has an automatic imbalance shut off system.

You can only start the unit when the centrifuge lid is closed.

In case of operation malfunction, an error code will appear on the display; this helps to indicate the reason for the malfunction and to eliminate it.

## Efficient and Versatile

The small, compact design is perfect for tight spaces.

Safety of your samples: even after a long centrifugation time, the temperature of your samples is at most 6 Kelvin higher than the ambient temperature.

Even with the high speed output of 6,800 rpm; the working environment is vibration-free, with minimal noise output, due to its' high stability.

Sleep Mode Function: if the centrifuge lid is open and the unit is not in use for more than five minutes, the unit automatically goes into sleep mode (This function can be de-activated!).

## Z 207 A



CHF 1120.-  
incl. Rotor

+2,000 „Falcons  
15ml“ for free\*  
value 200.-



### technical specification

max. Speed	6800 rpm
max. RCF	4445 xg
max. Volume	8 x 15 ml RB / 4 x 15 ml Falcon
Speed range	200 - 6800 rpm
Running time	59 min 50 s / 10 s increments 99 h 59 min / 1 min increments
Dimensions	(WxHxD): 28 x 24 x 35 cm
Weight	8 kg

## Z 207 M



CHF 1500.-  
incl. Rotor

up to 10,000  
„Eppis“ for free\*  
value 300.-



### technical specification

max. Speed	13,500 rpm
max. RCF	17,317 xg
max. Volume	24 x 1.5/2.0 ml
Speed range	200 – 13,500 rpm
Running time	59 min 50 s / 10 s increments 99 h 59 min / 1 min increments
Dimensions	(WxHxD): 28 cm x 24 x 35 cm
Weight	12.5 kg

## Accessories for Z 207 M



Order No.	Angle Rotor	Angle
221.59 V01	24 x 1.5/2.0 ml	45°
221.59 V02	24 x 1.5/2.0 ml	45° hermetically sealable
221.58 V01	18 x 2.0 ml (Spin Column Kits)	45°

## \*Choose your Tubes

Z 207 A	
53015	Centrifuge Tubes 15ml PP, with HDPE cap, conical, graduated „Falcons“
Z 207 M	
675015S	Micro Tube 1,5ml with flat cap, transparent, Safe lock, (5000 for free) „Eppis“
61015	Micro Tube 1,5ml with flat cap, 5 colours assorted blue, pink, yellow, orange and green „Eppis“
60015	Micro Tube 1,5ml with flat cap, transparent cristal clear, Premium version „Eppis“
60115	Micro Tube 1,5ml with flat cap, steril, transparent, clear Premium version „Eppis“